


Capre MC Medical Chair

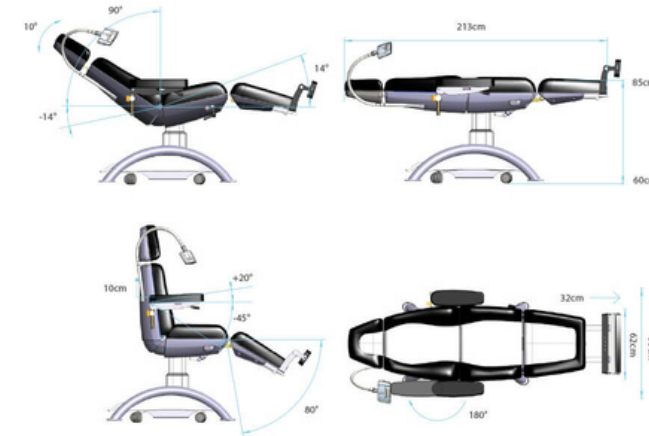

MC001 Medical Chair

The Capre MC Medical Chair is a premium, multifunctional treatment chair engineered with clinical precision for demanding healthcare settings such as dialysis, cytostatic therapy, post-day surgery recovery, and general examinations.

Collaboratively designed with healthcare professionals to adapt to evolving care environments, it offers a contemporary, ergonomic alternative to traditional hospital beds.

The chair includes pre-programmed positions, plus a user memory setting and a one-touch resuscitation option. Its Trendelenburg tilt and motorized backrest, seat, height, and leg adjustments are all controlled via a seamless handset or optional foot pedals. Mechanical adjustments extend the head and leg sections, further enhancing versatility.

With a contoured design, the Capre MC delivers high comfort via a insulation layer and meticulously shaped upholstery. Multi-function armrests remain level regardless of backrest angle, and can rotate 180°.

Code	Total Length	Height Adj	Total Width	Patient Surface Width	SWL
MC001-COLOUR	186-213cm	60-85cm	88cm	62cm	160kg

Code	Backrest Adj	Leg Section	Seat Tilt	Armrest Adj
MC001-COLOUR	-14° to +90°	-80° to 0°	+12° to 0°	up to 10 cm

Configuration

MC001-Colour: Capre MC Medical Chair, Inc USB Charger, Hand Control with Flexible Arm, Pre-programmed Positions & Motorised Leg Section



-FALLS



-ONYX



-TITANIUM



-BEACH

Accessories	Code
Plastic Cover for Leg Section	R19888
Plastic Cover for Foot Support	R19889
IV-Pole, Straight, One-Hand Height Adjustment, 4 Hooks	60121
Black Tray, Attached to the Upper Frame of The Chair	Coffee Tray
Flexible Tablet Holder	Tablet Support
Led Light / Lamp with Flexible Stand	60131L
Foot Control for Height & Back Section Adj	Foot control RC/MC
Fixed Foot Control for Trendelenburg Adj	Foot trend RC/MC
Neck Pillow	RCNP
Motorised Leg Section with Foot Support	MC19630
Battery Back-Up	BA21RC/MC



Handling

The Capre MC Medical Chair is designed to deliver seamless functionality across a variety of clinical settings, from outpatient clinics to specialist treatment rooms.

Its smooth-rolling castors and compact frame provide easy mobility, allowing staff to reposition the chair quickly between bays or departments.

The lightweight yet sturdy construction, combined with intuitive handset controls, makes height, tilt, and backrest adjustments effortless, helping clinicians optimise workflow while maintaining precision in patient positioning.

Comfort

With a focus on patient well-being during longer treatments and examinations, the Capre MC incorporates ergonomically designed cushions manufactured from premium, high-resilience foam. Upholstered in antimicrobial, soft-touch fabric, the chair provides both hygiene assurance and a welcoming feel for patients.

The adjustable backrest, seat, and leg sections are designed to reduce pressure points, promote circulation, and support a natural sitting posture. The overall design encourages relaxation while ensuring clinicians have unobstructed access for efficient procedures.

Safety

The Capre MC combines modern design with robust safety features, ensuring peace of mind for both patients and practitioners. Its durable frame and stable lift mechanism provide secure positioning throughout adjustments.

Integrated emergency positioning functions, including Trendelenburg, are easily accessed through the handset to allow immediate response in critical situations. All materials and components conform to the highest infection control standards, fire safety regulations, and European medical device directives, supporting reliable, compliant performance in healthcare environments.