


# Lifting Hoists

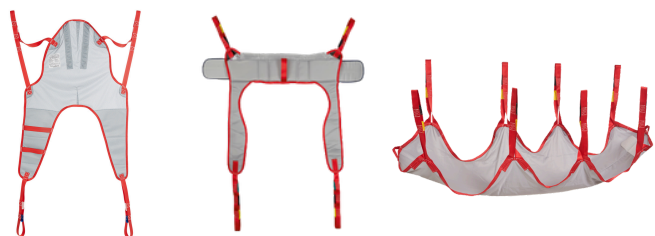

## Arnold & James Patient Hoists

Built with a high-strength epoxy-coated steel frame, the range delivers exceptional durability while remaining easy to manoeuvre thanks to its four swivel castors, two of which are braked for secure positioning. The low under-ride profile allows smooth access under most beds and furniture, while the base width is adjustable via a foot-operated mechanism to accommodate different patient situations and equipment layouts.

The Arnold 150 with Scale PW 600 combines all the core features of the standard model with an integrated medical-grade digital scale, enabling weight measurement during transfer to support accurate clinical assessment and care planning. All models are equipped with electric lifting powered by a rechargeable battery pack, emergency stop and mechanical lowering systems, ergonomic push handles, and compatibility with a range of slings, ensuring both patient comfort and operator safety.

Specifically designed to assist patients with rising from seated to standing positions, the James 150 Sit-to-Stand hoist is ideal for low mobility use cases. It shares the same epoxy-coated steel frame and portability via four swivel castors (two braked) as the Arnold, while offering tailored support through leg straps, 180° grab handles, and angled footrests.



Code	Lift Height Min.	Lift Height Max.	Underpin Height	Weight	SWL
310.00.10	55cm	174cm	14cm	43kg	125kg
320.00.10	55cm	174cm	14cm	43kg	150kg
320.06.10	55cm	174cm	14cm	43kg	150kg
320.02.10	84cm	150cm	14cm	53kg	150kg

Accessories	Code
Toilet Sling w/ Velcro Fastener Size M, L or XL	321.50.20/ 321.50.25/ 321.50.30
Net/Bath Sling Size S, M, XL/XXL	321.20.10/ 321.20.21/ 321.20.30
Back Sling w/ Velcro Fastener Size M/L, XL, XXL (James Only)	321.70.20/ 321.70.30/ 321.70.40
Wide Hanger- Four Point Spreader Bar	420.10.42
Horizontal Sling	321.60.30

## Configuration

- 310.00.10: Arnold 125 Patient Hoist
- 320.00.10: Arnold 150 Patient Hoist
- 320.06.10: Arnold 150 Patient Hoist w/ Scale PW 600
- 320.02.10: James 150 Standing Up Hoist

