

Capre RC Recovery chair

For recovery after outpatient surgery



Product pictured with optional accessories.



Watch a video
<http://bit.ly/capretuoli>



Main features:

1. Modern solution for post-day surgery recovery for improving vitality during recovery process
2. Comfortable to sit, lay down or relax in semi-upright position
3. Easy to use adjustments such as pre-programmed positions as well as a user memory position

The luxurious Capre multipurpose chair has been designed together with surgeons to fulfil the demands of today's healthcare processes in post-day surgery recovery. The chair is used instead of the traditional hospital bed to improve patient comfort during recovery. The Capre chair features pleasantly contoured upholstery, individual adjustment possibilities (with up to five motors) and many optional built-in comforts, such as a tablet pc stand and reading light. All of this combines to improve patients' quality of life during recovery, while also improving the quality image of the healthcare facility.

The control system allows the independent adjustment of the chair height, back section, seat angle and leg sections. The leg section is split and motorised. Both the head section and the split leg section also feature a length adjustment mechanism. The chair features a Trendelenburg position as well as three pre-programmed treatment positions (chair, beach chair and flat lying surface) and a user memory position in which the customer's own settings can be stored. If needed, the chair can also be set to either lying down or resuscitation position by one press of a button. Due to its modular design, the chair is also suitable as an examination chair.

LOJER® For easy care


Technical data and features	
Total length	186 - 213 cm
Patient surface length	186 - 213 cm
Total width	880 mm
Patient surface width	620 mm
Total weight	100 kg
Height adjustment	60 - 85 cm
Head section adjustment	10°, length adjustment 10 cm
Back section adjustment	-14°- +90
Back section length	600 mm
Seat adjustment	14°
Seat width	620 mm
Leg section adjustment	-80° - 0°
Arm support length	42 cm
Arm support adjustment	Height adjustment 10 cm, angle adjustment -45° to +20° and rotation 0° - 180°
Electric Trendelenburg	max. 14°
Central locking castors	Ø75mm
Other standard features	Split leg section with length adjustment USB charger Hand control with flexible arm or hand control with spiral cable Pre-programmed positions and one for user to decide White or grey cover for the base frame, under the seat and back section
Upholstery	Anatomically contoured Future -upholstery, several colour options
Safe Working Load (SWL)	160 kg
CE-marking	Yes
Warranty	2 years
Country of Manufacture	Finland

Factory installed options:

- Electric foot control for height and back section adjustment
- Electric Trendelenburg; fixed foot control for Trendelenburg adjustments
- All foot controls (featuring both foot control for height and back section + fixed control for Trendelenburg)
- Battery for back-up

Other accessories:

- Sample arm support with anatomic upholstery (length 60 cm)
- Neck pillow
- Breathing hole
- Black tray
- Flexible tablet holder
- LED light
- IV-pole

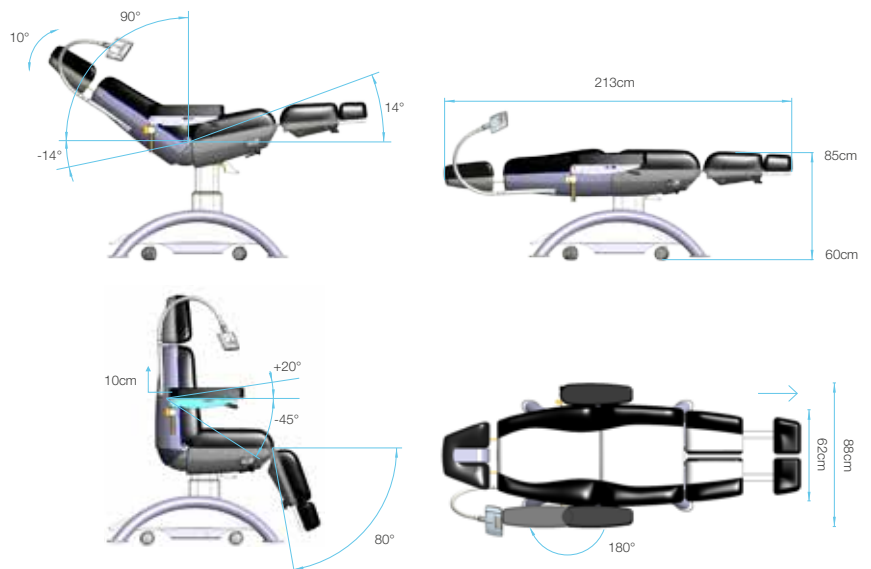

LOJER® For easy care